


**West** 

By your side  
for a healthier world™

WESTAR® RU AND RS STOPPERS AND PLUNGERS

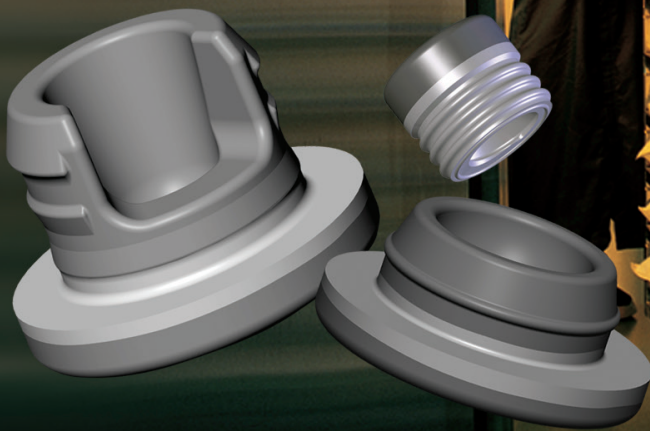

**Westar**®

Components Prepared for  
Pharmaceutical Use

# Flexible Solutions for Pharmaceutical Components



The West Ready Pack® system includes Westar ready-to-use stoppers.

Pharmaceutical manufacturers are challenged by ever-shortening production cycle times. A delay in having components prepared can ruin your manufacturing schedule.

To mitigate the risk of component preparation delays, you can apply Westar® RU (ready-to-use) or RS (ready-to-sterilize) components, including syringe plungers and stoppers. Westar components are washed with the validated Westar process.

## Choose Westar

Westar components can streamline your manufacturing operations by eliminating many of the processes related to component preparation. As a result, you may benefit from:

- Reduced testing and validation studies
- Increased throughput
- Reduced labor and maintenance costs
- Improved compliance
- A Drug Master File (DMF) to support chemistry, manufacturing and controls (CMC) filing

Westar components are packed in plastic cartons and may be placed on plastic shipping pallets to minimize particle contamination and facilitate handling within your controlled environment.

## Choose Ready-to-Use

For faster processing, choose RU components. Ready-to-use components are sterilized and delivered to you in integrity-tested packaging systems, which include rapid-transfer port bags, to ease introduction of components into your filling line.


With RU components, you minimize variability in your sterile operations and may even eliminate component preparation from your plant entirely.

## Choose Ready-to-Sterilize

When choosing RS components, you select clean components that comply with applicable USP and Ph. Eur. compendial standards and provide consistency of preparation from clinical trials through commercialization. West performs a portion of your component preparation process.

Our testing and validation work can help improve the speed of product development and manufacturing.

## Westar RU and RS Plungers and Stoppers

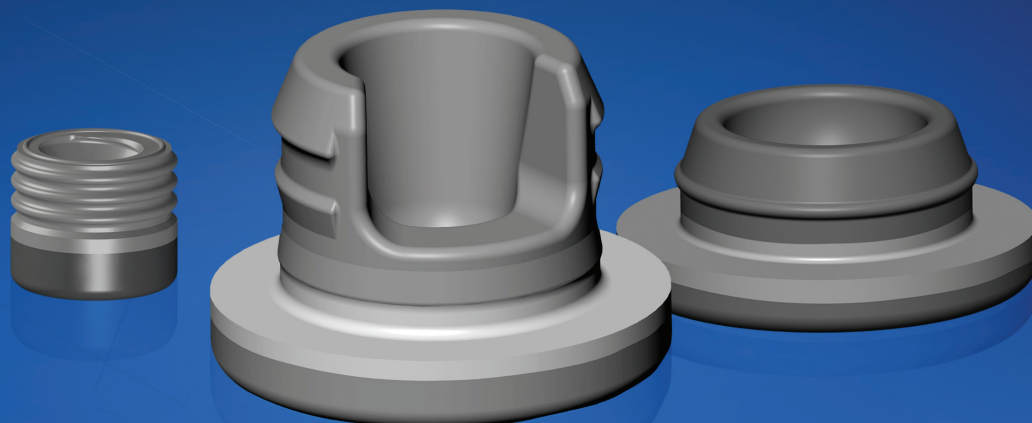

## Westar Component Processing

The validated process ensures that components are washed and rinsed with Water For Injection and packed in an ISO 5 clean room.

### Validated for Quality Assurance

Westar components are prepared using a documented, validated process that helps companies comply with current regulatory requirements. All product lots are tested and accompanied by a Certificate of Analysis as documented evidence of compliance to specifications, such as:

- Endotoxin
- Bioburden
- Particulate
- Silicone level (if applicable)

### Essential Compliance

The process complies with current FDA thinking on cGMP lot-to-lot closure preparation and testing. The wash process is supported by a DMF filing in the U.S. and Canada.

Since 1995, the Westar RS DMFs have successfully supported our customers' applications. To request a Letter of Authorization to reference the DMF in support of regulatory filings, use the on-line form at [westpharma.com/services/loa\\_form.asp](http://westpharma.com/services/loa_form.asp) or contact your West representative.

### Risk Mitigation

Westar components are available from West's facilities in North America, Europe and Singapore. The facilities adhere to applicable cGMP requirements and produce components with specified bacterial endotoxin, bioburden and particulate levels for primary pharmaceutical closures.

### The Leading Choice, The World Over

The world's leading pharmaceutical companies choose Westar components to help ensure consistency of preparation from clinical laboratory testing through commercial production. Call your West

representative to learn how Westar components can be easily specified into your operations.

### Select a Complete Packaging System

For assurance of component fit and optimal container closure integrity, select the West Ready Pack® system. The system includes West's RU Flip-Off® seals that have been sterilized by a validated gamma irradiation or steam autoclave process. The seals are delivered in packages designed to fit into sterile filling and capping lines. Ready-to-use vials and stoppers in small quantity packaging complete the Ready Pack system.

West's Flip-Off® CCS – clean, certified, sterilized – seals, designed for aseptic operations, can help achieve compliance with the European Commission Guide to GMP, Annex 1.



By your side  
for a healthier world™

Every day, injectable drugs are administered to improve the lives of millions of patients around the world. And every day, West is working by your side to design and manufacture drug packaging and delivery systems that will bring your drugs from concept to the patient more efficiently, reliably and safely. West understands your challenges and helps with solutions every step of the way, with cutting-edge production technologies, an unmatched expertise in global regulatory compliance, and an ever-growing knowledge base of pharmaceutical drug product testing, development, packaging and delivery. Whether you want an end-to-end solution or your focus is on one piece of the process, West is by your side for a healthier world.

To learn how components with Westar® RU and RS processing can help reduce the risk of product loss and failures, call your West account manager or a West technical representative. We'll help you determine the solution that best meets your needs.

North America +1 800 231 3000

South America +55 11 4055 606

Europe +49 2403 7960

Asia Pacific +65 6862 3400

[www.WestPharma.com](http://www.WestPharma.com)

West and the diamond logo, Flip-Off®, Ready Pack® and Westar® are registered trademarks of West Pharmaceutical Services, Inc. in the United States and other jurisdictions. By your side for a healthier world™ is a trademark of West Pharmaceutical Services, Inc. in the United States and other jurisdictions. Copyright © 2013 West Pharmaceutical Services, Inc.

#7252 • 0413